# **IHDE**DENTAL ?





Dr. Ihde Dental has been a reliable partner for over 60 years providing a wide range of implant systems and consumables. We supply dentists and dental technicians with precisely coordinated materials and systems, which are easy and reliable to use. We always ensure high quality and an excellent price-performance ratio so that you can guarantee allround treatment for your patients that is cost-effective and highly efficient. The following catalog gives you an overview and all the essential information about our implant systems. You can also contact us personally any time using the phone numbers provided. Further information can be found on our websites:

## www.implant.com | www.ihde-dental.de | www.ihde.com

**The company** was founded in 1954 in Berlin by the dental technician Klaus Ihde. The company relocated to Bavaria in the 1960s. At the end of the 1980s, Dr. Ihde Dental GmbH (Germany) and Dr. Ihde Dental AG (Switzerland) were formed from the Klaus Ihde retail company. Ihde Dental is now represented in four locations in Europe and over 45 countries. The company group is one of the most innovative implant companies in the world – based on new developments and patents issued or pending.

**The core activities** of Ihde Dental are the development, procurement and distribution of medical products. We use a large number of suppliers in consumables, but we have produced implants in our own factory for many years. All components are manufactured quickly, precisely and economically thanks to state-of-the-art production technology and well-equipped machinery.

#### **Our partners**

Users and customers provide us with many new ideas and excellent suggestions. Collaboration with our customers is extremely important to us. Contact us at any time if you have any improvements or questions. Your ideas and opinions help us all to meet the daily wishes of patients to a greater and better extent. We also put the needs of the patient first..

#### Our market performance and work ethic

Since it was founded, the company has focused on innovative ideas and advanced technology, premium quality, an excellent price-performance ratio, optimal patient and user friendly products and durability. Our range combines the latest findings from research and practices in many countries around the world.

## Customer orientated to us means – available for you!

- We provide training courses, refresher courses and user advice.
- We provide customers with comprehensive and technically sound advice.
- We also visit you in your practice upon request.

Please call us to arrange an appointment or send us an email.



## THE ADVANTAGES

OF THE ENDOSSEOUS DENTAL IMPLANT SYSTEM BLP®

The surface of **Bone Level Plus®** implants provide a specially lasered surface with exactly defined properties. For anti-rotation an internal square connects with press-fit to the abutment. The cone in combination with the internal stare provides stability and 100% tightness. **Bone Level Plus®** implants are universally applicable for fixed and removable prosthetics.

The prescribed or recommended tightening torques for implants, abutments and screws can be found on our website:

www.implant.com/en/downloads

/ Safely ``antirotational thanks to its internal ``precision square /

Cone technology for a tight seal


Universally suitable for fixed and removable prosthodontics

Made of highly resistant titanium alloy

Smart instrument tray



BLP®

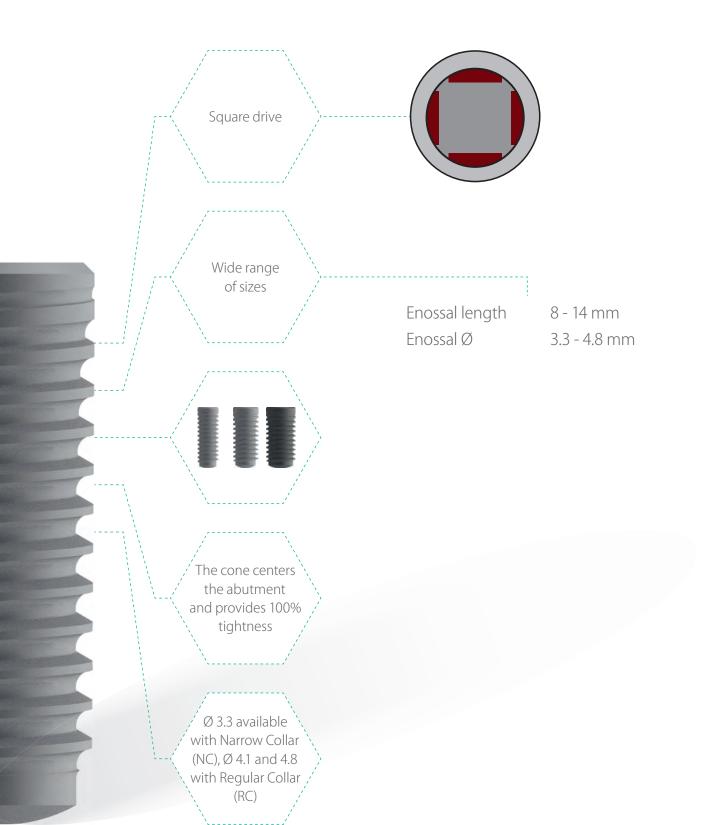


#### PREPARATORY STEPS WHEN USING A DRILLING TEMPLATE

- 1. Have your laboratory produce a drilling template with the appropriate drill holes for the marker bore. To be on the safe side, the laboratory might insert guide sleeves (**REF** BFH) into the drill holes to ensure that the drilling angle is correct. Use a 2 mm ø drill for pilot drilling.
- 2. For subsequent drilling sequences, drill stops can be used that are slid over the drill according to the appropriate depth of the drill hole and screwed in place. Consider the thickness of the mucosa and the height of the template as appropriate.

Thanks to the extremely high cutting efficiency of our drills, no ascending drilling sequences will usually be required.


Recommended RPM: 2000-5000

Apply sufficient cooling and allow the cooling to reach the working blades of the drills. Drill stop taking from  $\emptyset$  2.8.

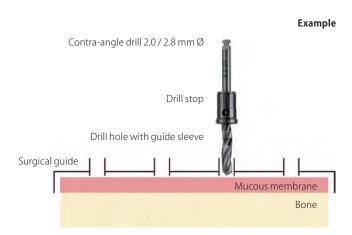

## **SURGERY**

## 1. Recommended drill sequence

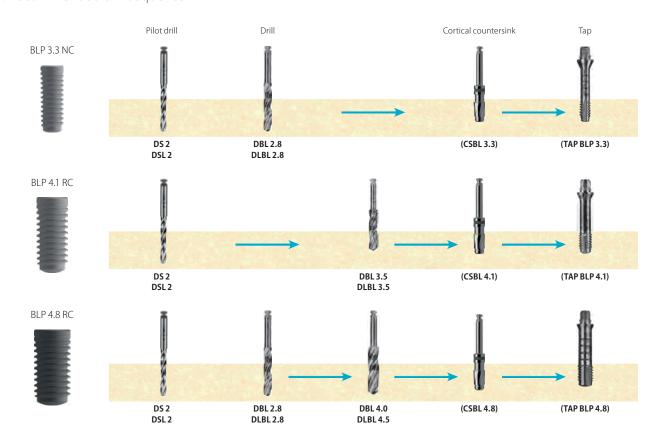

Owing to the high quality and geometry of the blades of our drills, the final preparation may be performed immediately after the pilot drilling.

BLP®

## 2. Implant packaging



Original packaging



Open the sealed cover at the lid. Remove the label and place it into the patients record.



The open pack contains the implant in a sterile tube (primary packaging).

## 3. Remove the implant from its packaging

- 1. Open the Lid
- 2. The implant is attached to the cap and can be removed by breaking it off at the pre-determined breaking point
- 3. Remove the implant, making sure not to touch the inner wall of the tube




## 4. Handling


#### 4.1 Connect

Attach the placement aid to the implant, holding the cap to which the implant is attached with the other hand.

#### 4.2 Mounting the adapter ITV WST / contra-angle

Place the ITV Wst (angled handpiece) or IT ITV (ratchet) adapter on the ITV BLP placement aid. Mount the placement aid. Hold the cap firmly in one hand and break off the implant at the predetermined breaking line.

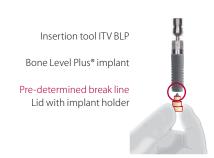

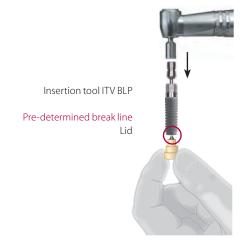

#### 4.3 Alternative to 4.2:

Place the IT ITV (ratchet) adapter on the ITV BLP placement aid.

Mount the placement aid. Hold the cap firmly in one hand and break off the implant at the pre-determined breaking line.



#### **5.** Insertion

Use the angled handpiece, ratchet or shank to screw the implant into the implant bed (clockwise).

The enossal aspect of the implant must be submerged in the bone. Upon **complete** insertion, the implant may be turned back ¼ revolution to reduce the load on the bone.

The system is suitable for deep insertion (below bone level).



## 6. Remove insertion tool from implant

Loosen the insertion aid from the implant by pulling it off.




## 7. Result




## 8. Aftercare

Seal the implant with a matching cover screw.



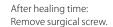

Arrow: Top view. Rotation counter clockwise.

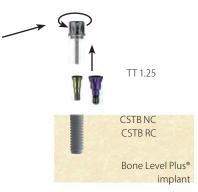

BLP®

9

## **9.** Pick-up impressions

#### **9.1** Impression with perforated custom tray



Insert impression posts HLT BLP NC/RC

Bone Level Plus® implant




#### **9.2** Before taking the impression

Take an impression in an A silicone. You can use the open-tray or the closed-tray technique.

It is necessary to remove the HLT BLP NC/RC impression post from the implant to be able to take out the impression tray.



Bone Level Plus® implant

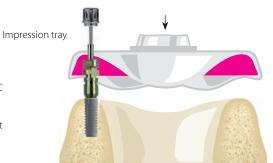


#### **9.3** Taking the impression

Detach the HLT BLP NC/RC from the implant. HLT BLP NC/RC will stay within the impression. Use TT 1.25 to loosen screw

Relief window in the impression tray


HLT BLP NC/RC

Bone Level Plus® implant

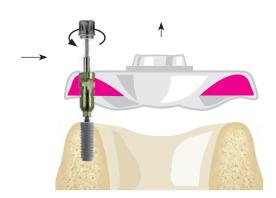

#### 9.4

View of the impression post in the impression (pick-up technique, bottom view).

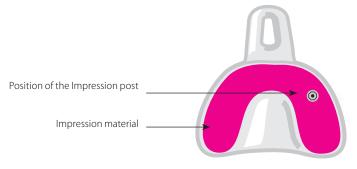

#### 9.5

Once the impression has been taken, the implant is closed with a healing abutment, while the impression is sent to the laboratory.

TT 1.25

Arrow: Top view. Rotation clockwise

Surgical screw CSTB NC/RC

Bone Level Plus\* implant

## 10. Closed tray impression taking

#### **10.1** Impression with closed tray

Impression with custom tray.

Securing the impression post with the thumbscrew

TS BLP NC/RC

Bone Level Plus® implant




#### **10.2** Before taking the impression

Take an impression in an A silicone. You can use the open-tray or the closed-tray technique.

With the closed impression technique, the TS BLP NC/RC will always remain on the implant when removing the impression.

Impression post
TS BLP NC/RC

Bone Level Plus® implant

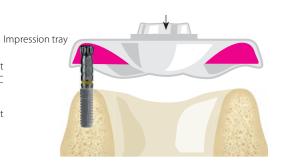


#### 10.3 Removing the impression

In the case of closed impressions, the TS BLP NC/RC impression post will remain on the implant after removing the impression tray.

The impression post will be removed afterwards.

Impression post TS BLP NC/RC

Bone Level Plus® implant

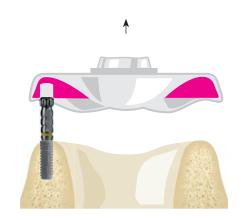

#### 10.4

Once the impression has been taken, the implant is closed with an HA NC/RC healing abutment, while the impression is sent to the laboratory.

TT 1.25 insert surgical screw CSTB NC/RC

Bone Level Plus® implant



## 11. Procedures in the laboratory


#### 11.1 Pick-up technique

Tighten the IAB against the HLT BLP (NC/RC) impression post.

Use the TT 1.25 to insert the lab analogue

HLT BLP NC/RC

IAB NC or IAB RC

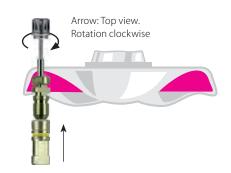

#### 11.2 Closed technique

Secure the IAB NC/RC against the TS BLP (NC or RC)  ${\color{red} (\rm A)}$ 

Reposition the impression post inside the impression  ${\color{red} {\bf 8}}$ 

Pour the impression.

Use the thumbscrew to tighten the impression post on the lab analogue.


TS BLP NC/RC

IAB NC or IAB RC



#### 11.3

Pour the impression in dental stone, then remove the impression posts from the lab analogues.

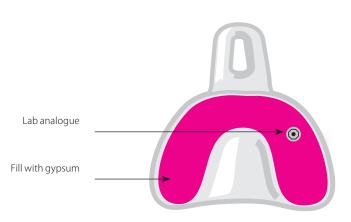

#### 11.4

The lab analogueue will now be embedded in the gypsum in the correct position.

IAB NC/RC



#### 11.5

Positioning of the screw-retained TLA2 15 BLP RC abutment, determining its optimal position and correct angulation.

**NOTE** The square end must be inserted completely into the analogue.

TT 1.25

Insert screw

TLA2 15 BLP NC/RC Watch out for the correct square end position

IAB NC/RC



#### 11.6

The correct position of the abutment must be ensured during transfer to the mouth.

TLA2 15 BLP NC/RC



#### 11.7

If multiple angled abutments are used, the laboratory will produce a removable resin splint (e.g. from pattern resin) to facilitate positioning within the mouth.

TLA2 15 BLP RC

Pattern Resin



## **BONE LEVEL PLUS® IMPLANTS**


With No-Itis® Laser Surface. The implant body is made of Ti6Al4V.









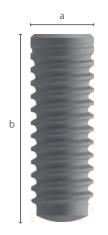

| a) Enossal Ø      | 3.3 - 4.8 mm  |
|-------------------|---------------|
| b) Enossal length | 8 - 14 mm     |
| NC                | Narrow Colla  |
| RC                | Regular Colla |

| Description   |  |
|---------------|--|
| BLP 3.3 8 NC  |  |
| BLP 3.3 10 NC |  |
| BLP 3.3 12 NC |  |
| BLP 3.3 14 NC |  |
| BLP 4.1 8 RC  |  |
| BLP 4.1 10 RC |  |
| BLP 4.1 12 RC |  |
| BLP 4.1 14 RC |  |
| BLP 4.8 8 RC  |  |
|               |  |


BLP 4.8 10 RC

BLP 4.8 12 RC

BLP 4.8 14 RC

| Enossal Ø | Enossal length | REF    | Price cat. |
|-----------|----------------|--------|------------|
| 3.3 mm    | 8 mm           | 900500 | Н          |
| 3.3 mm    | 10 mm          | 900501 | Н          |
| 3.3 mm    | 12 mm          | 900502 | Н          |
| 3.3 mm    | 14 mm          | 900503 | Н          |
| 4.1 mm    | 8 mm           | 900504 | Н          |
| 4.1 mm    | 10 mm          | 900505 | Н          |
| 4.1 mm    | 12 mm          | 900506 | Н          |
| 4.1 mm    | 14 mm          | 900507 | Н          |
| 4.8 mm    | 8 mm           | 900508 | Н          |
| 4.8 mm    | 10 mm          | 900509 | Н          |
| 4.8 mm    | 12 mm          | 900510 | Н          |
| 4.8 mm    | 14 mm          | 900511 | Н          |

Min. Insertion torque 35 Ncm

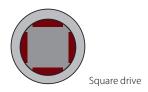

Delivery inclusive insertion tool ITV BLP and surgical screw REF 900518 or 900519





- Safely anti-rotational thanks to its internal precision square
- Cone technology for a tight seal
- Universally suitable for fixed and removable prosthodontics
- The cone centers the abutment and provides 100% tightness

## **SURGICAL SCREWS**



| Description                        | Code    | REF    | Price cat. |
|------------------------------------|---------|--------|------------|
| Surgical screw for BLP 3.3         | CSTB NC | 900518 | В          |
| Surgical screw for BLP 4.1 and 4.8 | CSTB RC | 900519 | В          |

## **GINGIVAFORMER**







| Description   |              | Code           | REF    | Price cat. |
|---------------|--------------|----------------|--------|------------|
| Gingivaformer | conical      | GF NC 3.6 2    | 900590 | В          |
| Gingivaformer | conical      | GF NC 3.6 3.5  | 900591 | В          |
| Gingivaformer | conical      | GF NC 4.8 3.5  | 900594 | В          |
| Gingivaformer | conical      | GF RC 4.5 2    | 900596 | В          |
| Gingivaformer | conical      | GF RC 4.5 4    | 900597 | В          |
| Gingivaformer | conical      | GF RC 4.5 6    | 900598 | В          |
| Gingivaformer | conical      | GF RC 62       | 900599 | В          |
| Gingivaformer | bottle shape | GFB NC 3.3 3.5 | 900602 | В          |
| Gingivaformer | bottle shape | GFB NC 3.3 5   | 900603 | В          |
| Gingivaformer | bottle shape | GFB RC 4.4 4   | 900604 | В          |
| Gingivaformer | bottle shape | GFB RC 4.7 6   | 900605 | В          |

## **BUR CYLINDER**



**Description**Bur cylinder for BLP 3.3 for telescope crowns
Bur cylinder for BLP 4.1 and 4.8 for telescope crowns

Recommended insertion torque  ${\bf 30\,Ncm}$ 

 Code
 REF
 Price cat.

 FZB NC
 900524
 D

 FZB RC
 900527
 D

## **ANALOGUES**





| Description                          | Code      | REF    | Price cat. |
|--------------------------------------|-----------|--------|------------|
| Implant analogue for BLP 3.3         | IA BLP NC | 900525 | В          |
| Implant analogue for BLP 4.1 and 4.8 | IA BLP RC | 900526 | В          |

15

## **STANDARD ABUTMENTS**

|   | Description   | Code     | REF    | Price cat. |
|---|---|----------|--------|------------|
|   | Abutment for cementing on BLP 3.3, step 1 mm high Height above step 4 mm, incl. screw SFBC NC           | CAB 1 NC | 900554 | E          |
| 1 | Abutment for cementing on BLP 3.3, step 3 mm high Height above step 4 mm, incl. screw SFBC NC           | CAB 3 NC | 900555 | E          |
| ₩ | Abutment for cementing on BLP 4.1 and 4.8, step 1 mm high Height above step 5.5 mm, incl. screw SFBC RC | CAB 1 RC | 900551 | E          |
| 1 | Abutment for cementing on BLP 4.1 and 4.8, step 3 mm high Height above step 5.5 mm, incl. screw SFBC RC | CAB 3 RC | 900552 | Е          |
| Ψ | Recommended insertion torque 20 Ncm   |          |        |            |

## **SCREW-RETAINED ABUTMENTS** (REDUCIBLE, GRINDABLE)



| Description  | Code           | REF    | Price cat. |
|--|----------------|--------|------------|
| Abutment Incl. screw SF B  | TAB BLP NC/RC  | 900521 | D          |
| Abutment for BLP 3.3, 15° angled<br>Anti-rotational, incl. screw SFB NC      | TLA2 15 BLP NC | 900528 | F          |
| Abutment for BLP 4.1 and 4.8, 15° angled Anti-rotational, incl. screw SFB RC | TLA2 15 BLP RC | 900523 | F          |
| Recommended insertion torque 20 Ncm  |                |        |            |

## **ANATOMICAL ABUTMENTS**



| Description  | Code    | REF    | Price cat. |
|--|---------|--------|------------|
| Anatomical abutment for BLP 3.3<br>Anti-rotational, incl. screw SFB NC         | ANAB NC | 900544 | F          |
| Anatomical abutment for BLP 4.1 and 4.8<br>Anti-rotational, incl. screw SFB RC | ANAB RC | 900543 | F          |

Recommended insertion torque  ${\bf 20~Ncm}$ 

## TITANIUM BASE FOR CAD CAM



| Description   | Code      | REF    | Price cat. |
|---|-----------|--------|------------|
| Titanium base for BLP 3.3, anti-rotation Incl. screw SFB NC | MB BLP NC | 900560 | D          |
| Titanium base for BLP 4.1 and 4.8, anti-rotation            | MB BLP RC | 900562 | D          |

## **CASTABLE ABUTMENTS**



| Description  | Material       | Code         | REF    | Price cat. |
|--|----------------|--------------|--------|------------|
| Castable abutment for BLP 3.3<br>Incl. metal base and screw      | CoCrMo/plastic | PLAB2 BLP NC | 900621 | G          |
| Castable abutment for BLP 4.1 and 4.8 Incl. metal base and screw | CoCrMo/plastic | PLAB2 BLP RC | 900623 | G          |

## PICK-UP IMPRESSION POST FOR PICK-UP IMPRESSIONS



| <b>Description</b> Impression post for BLP 3.3 | Code<br>HLT BLP NC | <b>REF</b> 900584 | Price cat. |
|--|--------------------|-------------------|------------|
| Impression post for BLP 4.1 and 4.8            | HLT BLP RC         | 900585            | С          |

## **IMPRESSION POST** FOR CONVENTIONAL IMPRESSIONS



| <b>Description</b> Impression post for BLP 3.3 | Code<br>TS BLP NC | REF<br>900586 | Price cat. |
|--|-------------------|---------------|------------|
| Impression post for BLP 4.1 and 4.8            | TS BLP RC         | 900587        | С          |
| Impression post long for BLP 3.3               | TSL BLP NC        | 900588        | С          |
| Impression post long for BLP 4.1 and 4.8       | TSL BLP RC        | 900589        | С          |

#### **ABUTMENTS FOR SCREW-ON PROSTHETIC**



| <b>Description</b> Gingiva height 0.5 mm | Code<br>TCT BLP NC 0.5 | REF<br>900635 | Price cat. |
|--|------------------------|---------------|------------|
| Gingiva height 1.5 mm                    | TCT BLP NC 1.5         | 900636        | D          |
| Gingiva height 3.5 mm                    | TCT BLP NC 3.5         | 900637        | D          |
| Gingiva height 0.5 mm                    | TCT BLP RC 0.5         | 900632        | D          |
| Gingiva height 1.5 mm                    | TCT BLP RC 1.5         | 900633        | D          |
| Gingiva height 3.5 mm                    | TCT BLP RC 3.5         | 900634        | D          |
|  |                        |               |            |

## Tighten with **HT 1.77**

## **IMPRESSION TAKING AND LABORATORY ACCESSORIES**

In this approach the position of the TCT hex is assigned.



|            | Transfer post | Long screw Tighten with HT 1.25 | TCT analogue | Castable<br>abutment<br>12mm high<br>Internally round<br>Pack of 5 | Castable<br>abutment<br>12mm high<br>Internally edged<br>Pack of 5 | Fastening<br>screw<br>Tighten with<br>HT 1.25 |
|------------|---------------|---------------------------------|--------------|--|--|---|
| Code       | TST           | SFL                             | ВТТ          | PSTR (grey)  | PSTA   | SF  |
| REF        | 418147        | 420428                          | 418100       | 418124   | 418123   | 418151  |
| Price cat. | В             | В                               | В            | В  | В  | В   |

## **LOCALICER®** FOR REMOVABLE PROSTHETIC

If LOC abutments are used in the upper jaw, we recommend to place at least six implants and to splint them through prosthetics in a stable manner. Tighten with **HT 1.77**.



| Description                    | Height | Code         | REF    | Price cat. |
|--------------------------------|--------|--------------|--------|------------|
| Localicer* for BLP 3.3         | 2 mm   | LOC BLP NC 2 | 900539 | D          |
| Localicer*for BLP 3.3          | 3 mm   | LOC BLP NC 3 | 900606 | D          |
| Localicer*for BLP 3.3          | 4 mm   | LOC BLP NC 4 | 900607 | D          |
| Localicer* for BLP 4.1 and 4.8 | 2 mm   | LOC BLP RC 2 | 900540 | D          |
| Localicer* for BLP 4.1 and 4.8 | 3 mm   | LOC BLP RC 3 | 900608 | D          |
| Localicer® for BLP 4.1 and 4.8 | 4 mm   | LOC BLP RC 4 | 900609 | D          |

## **ACCESSORIES FOR LOCALICER®**



| Description                                    | Code   | REF    | Price cat. |
|--|--------|--------|------------|
| Analogue + impression set                      | AA LOC | 462337 | C          |
| Set with 5 caps + 1 housing (EXTERNAL PRODUCT) | NCS    | 462338 | D          |

#### Pull-off force

Yellow 600 g, Pink 1.200 g, Transparent 1.800 g, Violet 2.700 g Black has no retention and is designed for temporary solutions for up to one month

## **MULTI-UNIT** ABUTMENTS

Insertion of the angled MU2 abutments with  $\mathbf{HT}$  1.25. Insertion of the straight MU2S abutments with  $\mathbf{HT}$  1.77



| <b>Description</b> Abutment 17° angled Incl. screw SFB RC                                       | <b>Material</b><br>Ti6Al4V | Code<br>MU2 17 BLP RC | <b>REF</b> 900640 | Price cat. |
|---|----------------------------|-----------------------|-------------------|------------|
| Abutment 35° angled<br>Incl. screw SFB RC   | Ti6Al4V                    | MU2 35 BLP RC         | 900641            | L          |
| Abutment straight<br>Gingiva height 0.5 mm  | Ti6Al4V                    | MU2S 0.5 BLP RC       | 900642            | G          |
| Abutment straight<br>Gingiva height 1.5 mm  | Ti6Al4V                    | MU2S 1.5 BLP RC       | 900643            | G          |
| Abutment straight<br>Gingiva height 2.5 mm  | Ti6Al4V                    | MU2S 2.5 BLP RC       | 900644            | G          |
| Gingivaformer incl. SF MU2<br>Height above abutment shoulder 6 mm                               | Ti6Al4V                    | GF MU 2               | 418286            | С          |
| Localicer® incl. SF MU2<br>Height above abutment shoulder 6.7 mm<br>Use with NCS Set REF 462338 | Ti6Al4V                    | MU 2                  | 418287            | С          |
| Prosthetic screw for MU2  | Ti6Al4V                    | SFB RC                | 900532            | В          |


## **ACCESSORIES** FOR MULTI-UNIT ABUTMENTS



## **INSTRUMENTS**

|           | Description                     | Code             | REF                | Price cat. |
|-----------|---------------------------------|------------------|--------------------|------------|
| 00 2.0    | Pilot drill short/long 2.0 mm Ø | DS 2 / DSL 2     | 425001 /<br>425002 | D          |
|           | Pilot drill short/long 2.8 mm Ø | DS 2.8 / DSL 2.8 | 425005 /<br>425006 | D          |
|           | Form drill short 2.8 mm Ø       | DBL 2.8          | 900570             | Е          |
|           | Form drill short 3.5 mm Ø       | DBL 3.5          | 900571             | Е          |
| 1033-1-10 | Form drill short 4.0 mm Ø       | DBL 4.0          | 900572             | E          |
|           | Cortical countersink 3.3        | CSBL 3.3         | 900576             | D          |
|           | Cortical countersink 4.1        | CSBL 4.1         | 900577             | D          |
| 0081 T.0  | Cortical countersink 4.8        | CSBL 4.8         | 900578             | D          |
|           | Тар                             | TAP BLP 3.3      | 900579             | D          |
|           | Тар                             | TAP BLP 4.1      | 900580             | D          |
|           | Тар                             | TAP BLP 4.8      | 900581             | D          |

## **GUIDE JACKET**

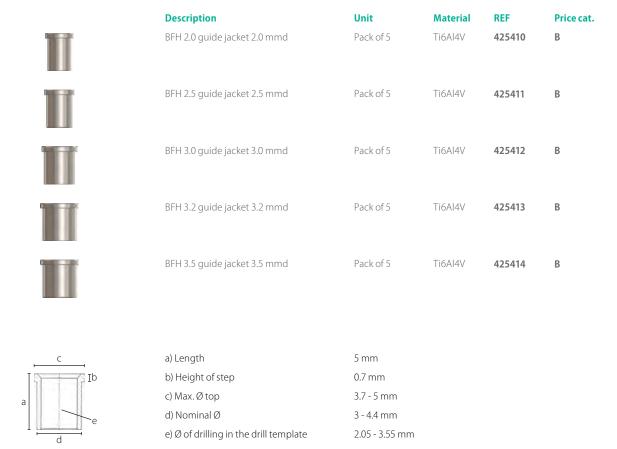

## **ADAPTER**



## **INSTRUMENTS** AND **TOOLS**

|        | Description   | Туре  | REF    | Price cat. |
|--------|---------------|---|--------|------------|
| O-OUT- | Ratchet RAT 2 | For all Hex instruments and insertion tools   | 425051 | K          |
|        | TW2           | Torque wrench, 10 - 70 Ncm. For all insertion tools, hex- and torxinstruments It is recommended to have the torque ratchets recalibrated by us once a year. | 425402 | S          |
|        | TT 1.25       | Torx instrument (for all screws)  | 425105 | С          |
| ====×  | TT 1.25 M     | Torx instrument (all screws) for contra-angle   | 425115 | С          |
|        | HT 1.77       | Hex instrument, long  | 425103 | С          |
|        | HTX 1.77      | Hex instrument, extralong   | 425104 | С          |
|        | PUW1          | Punch   | 425404 | С          |

## **STARTER** TRAY



| Description                 | Code        | REF      | <b>Price €</b> |
|-----------------------------|-------------|----------|----------------|
| Adapter contra-angle short  | ITVS        | 500851   |                |
| Adapter contra-angle medium | ITV M       | 500852   |                |
| Ratchet adapter for IT V    | IT ITV      | 500854   |                |
| Torx instrument             | TT 1.25     | 425105   |                |
| Cortical countersink 3.3    | CSBL 3.3    | 900576   |                |
| Cortical countersink 4.1    | CSBL 4.1    | 900577   |                |
| Cortical countersink 4.8    | CSBL 4.8    | 900578   |                |
| Pilot drill                 | DS 2.0      | 425001   |                |
| Form drill                  | DBL 2.8     | 900570   |                |
| Form drill                  | DBL 3.5     | 900571   |                |
| Form drill                  | DBL 4.0     | 900572   |                |
| Тар                         | TAP BLP 3.3 | 900579   |                |
| Тар                         | TAP BLP 4.1 | 900580   |                |
| Тар                         | TAP BLP 4.8 | 900581   |                |
| Torque wrench               | TW2         | 425402   |                |
| Starter tray empty          |             | 60045-K  | upon request   |
| Starter tray with content   |             | S60045-K | upon request   |

## **INSTRUMENT** TRAY



| Description                  | Code        | REF      | <b>Price €</b> |
|------------------------------|-------------|----------|----------------|
| Pilot drill                  | DS 2        | 425001   |                |
| Form drill                   | DBL 2.8     | 900570   |                |
| Form drill                   | DBL 3.5     | 900571   |                |
| Form drill                   | DBL 4.0     | 900572   |                |
| Standardized probe           | PDG         | 425400   |                |
| Standardized probe           | PDG         | 425400   |                |
| Standardized probe           | PDG         | 425400   |                |
| Cortical countersink 3.3     | CSBL 3.3    | 900576   |                |
| Cortical countersink 4.1     | CSBL 4.1    | 900577   |                |
| Cortical countersink 4.8     | CSBL 4.8    | 900578   |                |
| Тар                          | TAP BLP 3.3 | 900579   |                |
| Тар                          | TAP BLP 4.1 | 900580   |                |
| Тар                          | TAP BLP 4.8 | 900581   |                |
| Ratchet adapter for IT V     | IT ITV      | 500854   |                |
| Adapter contra-angle short   | ITV S       | 500851   |                |
| Adapter contra-angle medium  | ITV M       | 500853   |                |
| Adapter contra-angle long    | ITV L       | 500852   |                |
| Universal adapter            | UAW         | 425107   |                |
| Punch                        | PUW 1       | 425404   |                |
| Torx instrument              | TT 1.25     | 425105   |                |
| Drill extension              | DX 2        | 500704   |                |
| Torque wrench                | TW2         | 425402   |                |
| Instrument tray empty        |             | 60018-K  | upon request   |
| Instrument tray with content |             | S60018-K | upon request   |



(The products of this catalogue are CE marked (class I) and CE 1936 marked (class IIa and IIb) according to 93/42/EC Directive)

Commercial products that are not monitored by our notified body are declared as third-party products.

We are certified DIN EN ISO 13485, and annex II of EEC Directive 93/42 EWG (2007).

Product dimension described in this brochure may differ from reality for technical reasons.

Bone Level Plus® implants are protected by patents. Bone Level Plus® is a registered trademark.

In case that implants would be reprocessed (cleaned, resterilized) infections could occur, because no validated procedures for reprocessing are available.

#### Compilation and clarification of symbols on the pack:



Batch No.



Sterilized by gamma radiation



Non-sterile



Intended for use by dentists or surgeons only



Single use product



Instruction for use



Expiry date



Store in a dry place



Store tightly keep closed



Do not use if packing is damaged



Do not resterilize



Manufacturer



Production date



Catalogue number



Safely anti-rotational thanks to its internal precision square

Cone technology for a tight seal

Universally suitable for fixed and removable prosthodontics

The cone centers the abutment and provides 100% tightness

## **IHDE**DENTAL ?

#### Dr. Ihde **Dental AG**

Dorfplatz 11 CH - 8737 Gommiswald / SG Tel +41 (0)55 293 23 23 contact@implant.com www.implant.com

#### Dr. Ihde **Dental GmbH**

Erfurter Str. 19 D - 85386 Eching / Munich Tel +49 (0)89 319 761 0 Fax +49 (0)89 319 761 33 info@ihde-dental.de www.ihde-dental.de

